


Meas.Nr. 49026 **Date** Dienstag, 15 Mai 2012 02:34 **Personal Report No.:**
Material Yellow Paste 2 P6cl 80ml ZrO2 1mm 15min **M120102**
Description dedisp. in EtOH; 30s ultrasonic LABORETTE 17

Calculation Fraunhofer broad **TradeOff** 1.000000e+002
Scans Fine 100 **Scans Coarse** 0 **Channels** 51 **Beam Obscuration** 15.00 %
Meas. Range 0.08 [µm] - 42.30 [µm] **Pump** 6 **Ultrasonics** 0

%	< µm
5	0.11
10	0.15
15	0.20
20	0.26
25	0.31
30	0.37
40	0.47
50	0.59
60	0.72
70	0.87
80	1.08
90	1.41
95	1.71
98	2.08
99	2.33

%	< µm
4.0	0.1
14.6	0.2
42.6	0.5
51.2	0.6
65.6	0.8
76.4	1.0
84.4	1.2
91.9	1.5
95.9	1.8
97.6	2.0
99.5	2.5
99.9	3.0
100.0	3.5
100.0	4.0
100.0	5.0

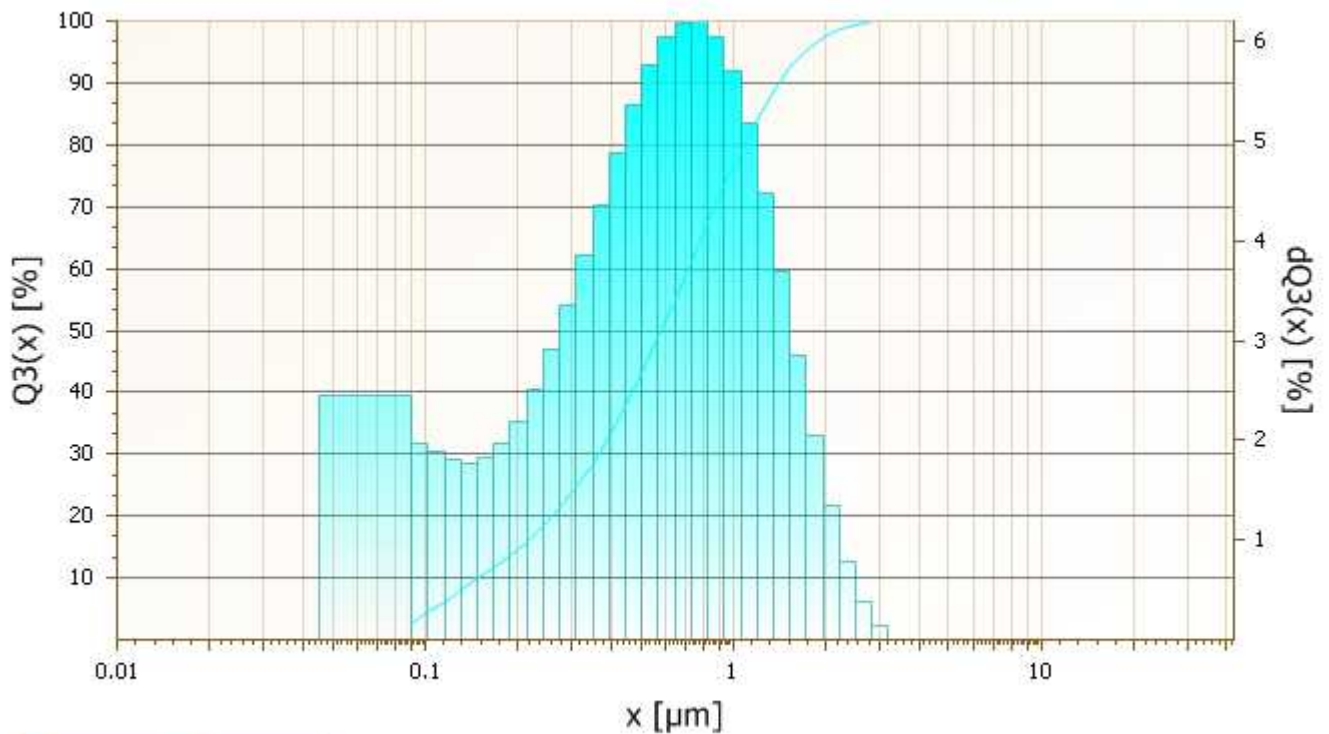

▲ 49026 dQ3(x) ▲ 49026 Q3(x)