


**Meas.Nr.** 62702      **Date** 04.04.2013 07:49      **Personal Report No.:**  
**Material** Zircon B1      120min P7pl 0,1mm ZrO2      **M130032**  
**Description** vordisp. in 0,1% Na4P2O7 mit 1min Ultraschall LABORETTE 17

**Calculation** Mie      n = 1.98970      -i 0      **broad TradeOff** 1.000000e+002  
**Scans Fine** 100      **Scans Coarse** 0      **Channels** 108      **Beam Obscuration** 13.00 %  
**Meas. Range** 0.01 [µm] - 42.30 [µm]      **Pump** 3      **Ultrasonics** 10

%	< µm
5	0.141
10	0.158
15	0.172
20	0.185
25	0.198
30	0.212
40	0.243
50	0.286
60	0.353
70	0.438
80	0.530
90	0.650
95	0.745
98	0.846
99	0.922

%	< µm
0.0	0.08
0.1	0.10
7.2	0.15
25.6	0.20
52.3	0.30
65.6	0.40
76.9	0.50
86.4	0.60
93.0	0.70
96.8	0.80
98.7	0.90
99.6	1.00
100.0	1.20
100.0	1.50
100.0	2.00

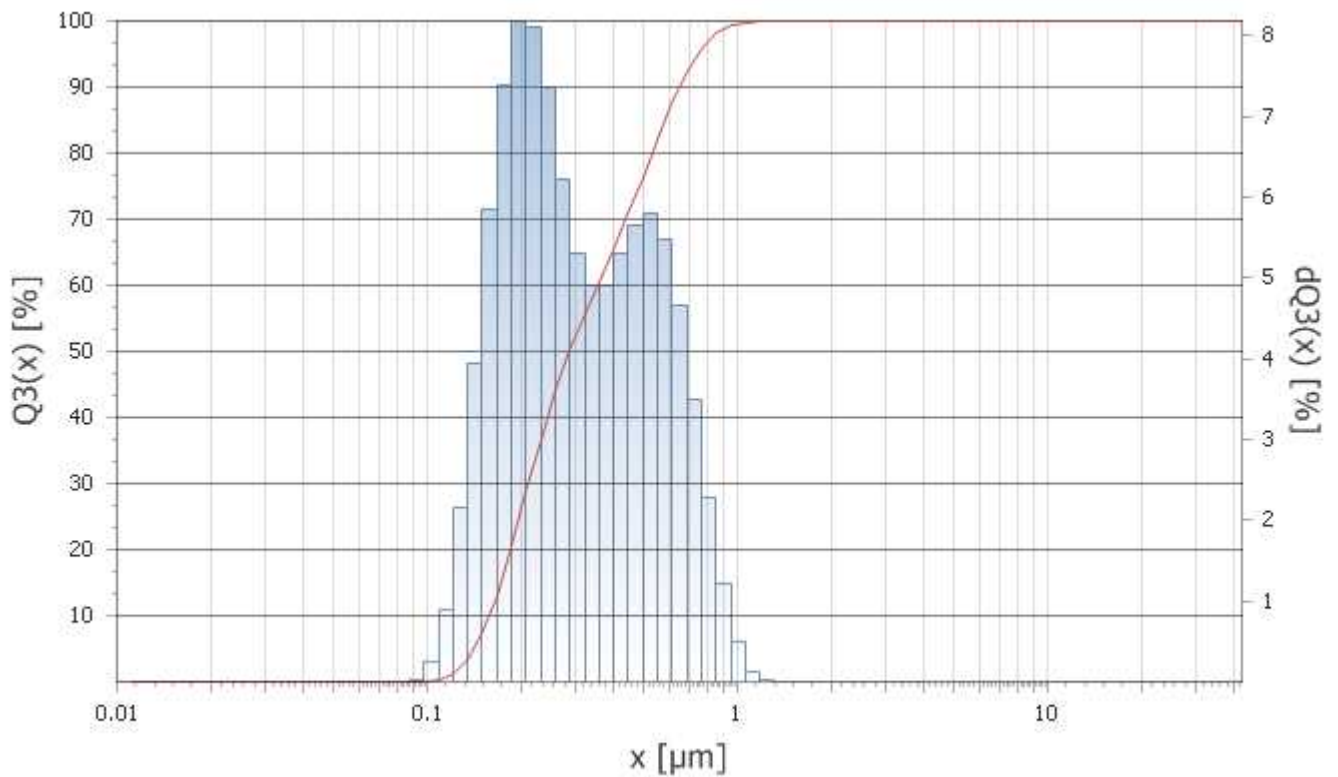

▲ 62702 dQ3(x)    — 62702 Q3(x)