


Meas.Nr. 62627 **Date** 04.02.2013 14:37 **Personal Report No.:**
Material Zircon B1 30min P7pl 2mm ZrO2 **M130032**
Description vordisp. in 0,1% Na4P2O7 mit 1min Ultraschall LABORETTE 17

Calculation Mie n = 1.98970 -i 0 **broad TradeOff** 1.000000e+002
Scans Fine 100 **Scans Coarse** 0 **Channels** 108 **Beam Obscuration** 16.00 %
Meas. Range 0.01 [µm] - 42.30 [µm] **Pump** 5 **Ultrasonics** 10

%	< µm
5	0.170
10	0.201
15	0.240
20	0.389
25	0.503
30	0.575
40	0.711
50	0.864
60	1.052
70	1.272
80	1.532
90	1.895
95	2.193
98	2.510
99	2.747

%	< µm
0.0	0.10
9.8	0.20
24.8	0.50
31.8	0.60
46.0	0.80
57.4	1.00
66.9	1.20
78.9	1.50
87.9	1.80
92.1	2.00
97.9	2.50
99.6	3.00
100.0	3.50
100.0	4.00
100.0	5.00

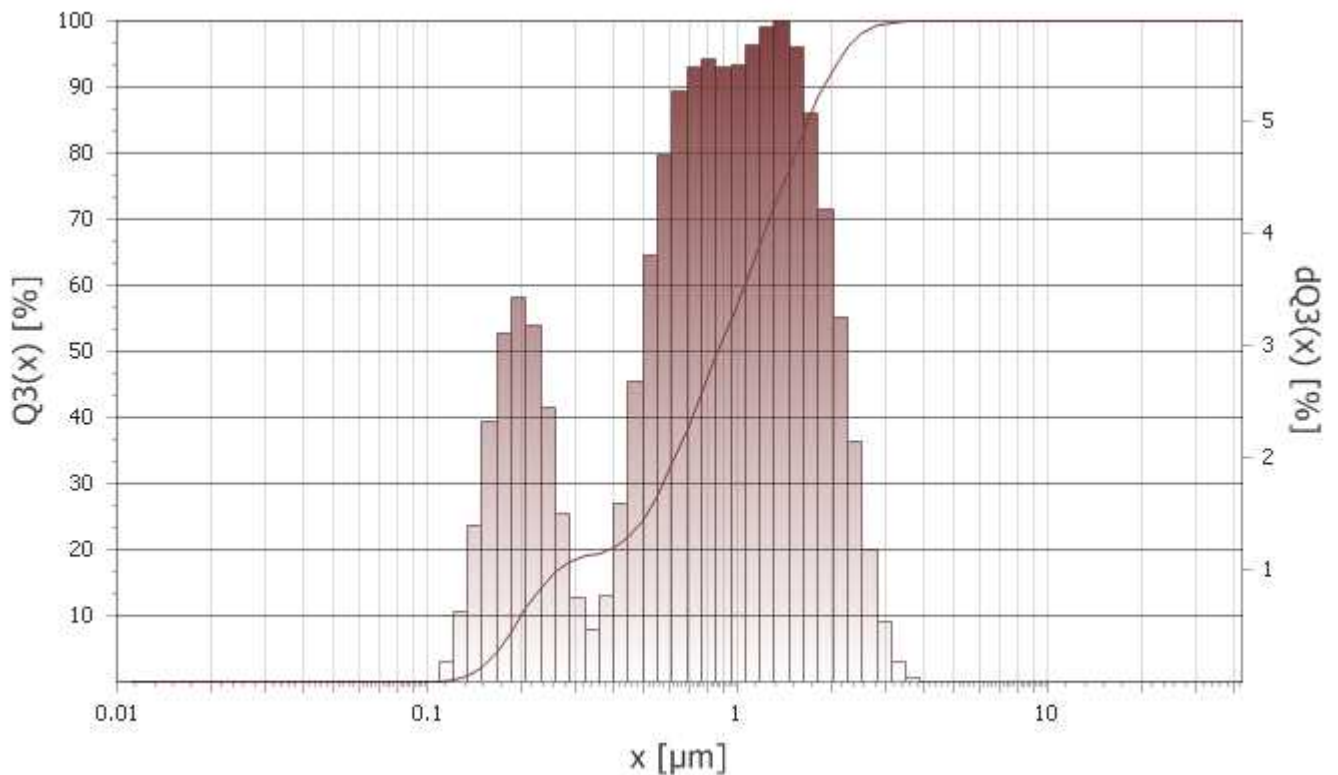

▲ 62627 dQ3(x) — 62627 Q3(x)