


Meas.Nr. 71091 **Date** 12.05.2013 13:38 **Personal Report No.:**
Material Asbestos 90 min P7pl 0,1 mm ZrO2 in water **M130287**
Description in water + 0,1 % Na4P2O7 - 1 min Ultrasonic LABORETTE 17

Calculation Mie n = 1.65 -i 0 **broad TradeOff** 1.000000e+002
Scans Fine 100 **Scans Coarse** 0 **Channels** 108 **Beam Obscuration** 11.00 %
Meas. Range 0.01 [µm] - 42.30 [µm] **Pump** 3 **Ultrasonics** 10

%	< µm
5	0.079
10	0.091
15	0.102
20	0.111
25	0.120
30	0.128
40	0.146
50	0.164
60	0.186
70	0.211
80	0.247
90	0.306
95	0.371
98	1.686
99	4.695

%	< µm
0.00	0.02
0.64	0.05
5.44	0.08
14.13	0.10
65.76	0.20
97.53	0.50
97.70	0.80
97.70	1.00
98.20	2.00
98.68	3.00
99.04	5.00
99.49	8.00
99.78	10.00
99.99	15.00
100.00	20.00

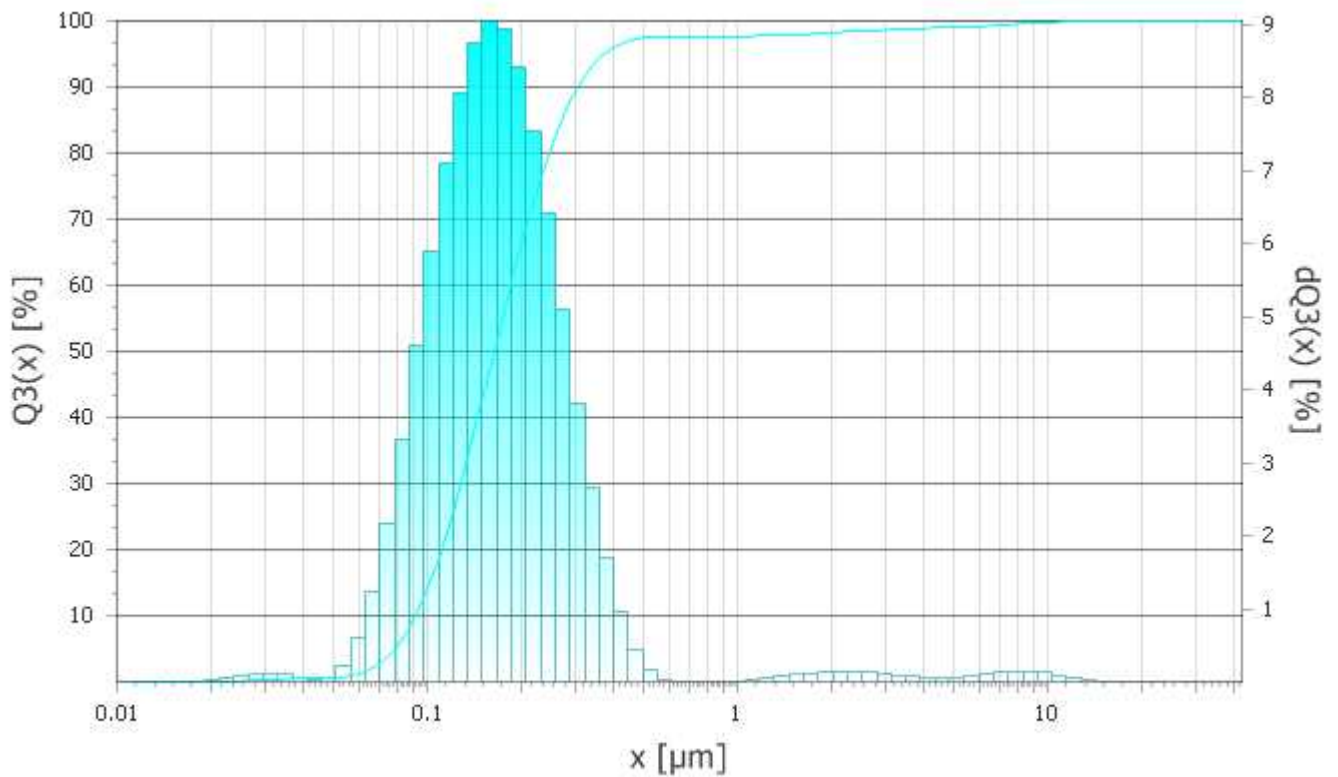

▲ 71091 dQ3(x) — 71091 Q3(x)